



HIGH FLOW BIO-X Filters

- air / gas filters
- PTFE impregnated borosilicate microfibre

HIGH FLOW BIO-X combines proven depth filter technology and a pleated construction to provide retention down to 0.01 micron in gas.

Flow rates typically 2-3 times that of membrane filters make HIGH FLOW BIO-X the filter that can dramatically reduce cartridge usage and installation size within the fermentation, food and beverage industries.

The specially developed PTFE impregnation process imparts greater strength and permanent hydrophobicity to the borosilicate microfibre media. This leads to excellent performance in applications such as the provision of sterile gas in filling machines.


Features and Benefits

- 94% voids volume PTFE impregnated microfibre
- Wide bore cartridge construction to maximize flow rate
- Stainless steel inner core
- Exceptionally high flow rates with low pressure drops
- Fully validated by aerosolized bacterial and viral challenge




Note: BIO-X is a registered trademark of Parker Hannifin Corporation.

Performance Characteristics

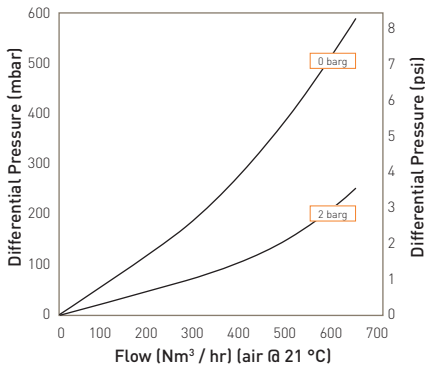

Flow rates for other sizes available upon request

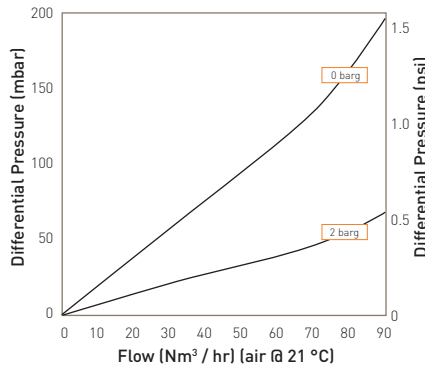

Flow rates for other sizes available upon request

10" Size (250 mm) Cartridge

A Size (125 mm) Cartridge

Specifications

Materials of Construction

- Filtration Media: PTFE Impregnated Borosilicate Microfibre
- Upstream Support: Polypropylene
- Downstream Support: Polypropylene
- Inner Support Core: 316L Stainless Steel
- Outer Protection Cage: Polypropylene
- End Caps: Encapsulated Polyethersulphone
- End Cap Insert: 316L Stainless Steel
- Standard o-rings/gaskets: Silicone

Biological Safety

Materials conform to the relevant requirements of 21CFR Part 177, EC1935 / 2004 and current USP Plastics Class VI - 121 °C and ISO10993 equivalents.

Recommended Operating Conditions

The maximum differential pressure in direction of flow (outside to in) is 3.0 barg (43.5 psig) at 70 °C (158 °F).

The maximum recommended continuous operating temperature is 70 °C (158 °F).

Effective Filtration Area (EFA)

10" (250 mm) 0.38 m² (4.09 ft²)

Sterilization

HIGH FLOW BIO-X cartridges can be in situ steam sterilized or autoclaved up to 142 °C (287.6 °F) for a maximum of 150 steam cycles.

Retention Characteristics

The HIGH FLOW BIO-X range of cartridges has been fully validated by aerosol bacterial challenge levels of 10¹² *Brevundimonas diminuta* per 10" (250 mm) filter cartridge. Independent test work also shows full retention to *MS-2 Coliphage*.

Integrity Test Data

All cartridges are integrity tested prior to despatch by the aerosol challenge test method using the Parker domnick hunter VALAIRDATA II.

Ordering Information

ZCHB - - -

Code Length (Nominal)	Code Endcap (10")	Code O-rings
B* 2.5" (65 mm)	C P-7	E EPDM
A* 5" (125 mm)	P BIO-X Retrofit	S* Silicone
K 5" (125 mm)	H UF Retrofit	V Viton
1 10" (250 mm)		
2 20" (500 mm)		
3 30" (750 mm)		
Code Endcap (Demi)		
H UF Retrofit		
T TRUESEAL		
Y Demi MCY		
Z Demi A & B Std		

* Supplied in packs of 3.

* Silicone o-ring supplied as standard without having to specify the 'S' code.



FILCO®, spol. s r.o.
 Dvorská 464/103, CZ-503 11 Hradec Králové, tel.: +420 495 436 233, info@filco.cz, www.filco.cz

Europe: ☎ +44 (0)191 4105121 📧 dhprocess@parker.com - North America: ☎ +1 608 824 0500 📧 dhpsales.na@parker.com 🌐 www.parker.com/dpharma

Parker domnick hunter has a continuous policy of product development and although the Company reserves the right to change specifications, it attempts to keep customers informed of any alterations. This publication is for general information only and customers are requested to contact our Process Filtration Sales Department for detailed information and advice on a products suitability for specific applications. All products are sold subject to the company's Standard conditions of sale.